

OPTICLUX 10



iNOVO[®]
AWARDS
2017



Details shown in the best light. OPTICLUX 10

The OPTICLUX 10 magnifier light combines a robust magnifier with outstanding light quality and ergonomic operation. As a certified medical product, it is used in dermatology, in clinics and hospitals for human, dental and veterinary medicine, in laboratories and veterinary medicine, in laboratories, but also for central steril supply department (CSSD).



OPTICLUX 10	technical data
luminous power *	6,000 Lux @ 0.15 m
colour temperature *	<ul style="list-style-type: none"> • 6500 K • optional: 365 nm (Wood's Light)
lens dimension	Ø 16 cm
colour rendering index	Ra > 95
usage	membrane switch on light head
class of protection	<ul style="list-style-type: none"> • version with cover: I • version with pin: II
protection rating	IP 20

special features	
OPTICLUX 10	<ul style="list-style-type: none"> • dimmable • maintenance-free • Wood's Light mode (optional) • lens cover • 1.9x magnification
variations	Further connector variants are available on request.
accessories	<ul style="list-style-type: none"> • table clamp • universal mount • fasting hardware • wall mount • roller stand

OPTICLUX 10 – variations



OPTICLUX 10 C
ceiling version



OPTICLUX 10 P
pin version



Product highlights



Optimally protected
sturdy plastic cover
protects the lens from
scratches and sunlight




360° turnover axis
for maximum light
availability



Especially for dermatology
two illumination modes in
one device:
optional Wood's Light for
examination of
fluorescent fungal or
bacterial infection on the
skin or scalp and
pigmentary changes.

Detailed description
to Wood's Light on page 110.

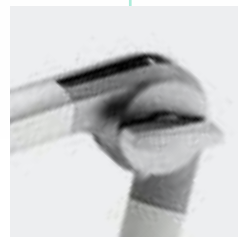

Ergonomic handling
spring arm and internal
friction joints allow
comfortable handling
without any risk of the
light head drifting



Flexible positioning
long lamp arm makes a
large radius of action up
to 142 cm possible